

# HA90C

Manual and Air-Drop Tables



  
Hill Laboratories  
COMPANY

# HA90C™

## MANUAL DROP TABLE

Designed with both you and your patient in mind, the standard HA90C electrically adjusts in height from 22" to 35" providing the ideal working height for you and easy access for your patients. The HA90C is also versatile—choose from the many available options to customize the HA90C for your specific needs.



*Shown in Gravel Upholstery with Beige Base.*

*Pictured with several options*


### ▼ STANDARD FEATURES

- Armrests, 11" face cut-out and paper roll
- Electrically-controlled height - 22" - 35"; Optional 18" to 32-1/2"
- Width - 24"; Length - 5'10" with drops; 6' without drops
- Foam Top - 2-1/2" Medium (standard); Soft or Firm
- Choice of 20 standard colors (or upgrade to more durable Finesse upholstery)
- Beige, Gray or Black Base, enclosed for safety
- One-year warranty on all parts and labor

### TOP PICK

The model shown is one of our most popular configurations with optional tilt, raised and dual-drop headpiece, thoracic drop, lumbar drop, pelvic drop and slide-out ankle support.

For even greater speed and convenience, consider the Air90C.

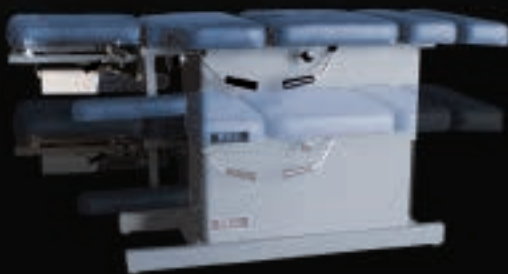

*The HA90C has electric adjustable height from 22" to 35" and available in three base colors; Beige, Gray and Black.*

# HA90C OPTIONS

Made-to-Order

## ▶ MANUAL OR AIR-DROPS


The standard HA90C can be equipped with up to 4 manual drops including cervical (straight-down and forward-motion), thoracic, lumbar and pelvic. The AirHA90C comes standard with all 4 air drops.

## ▶ THORACIC SPRING BREAKAWAY

The Thoracic Spring Breakaway (with release) allows for a thoracic recoil movement, features adjustable tension and locks at any point down to 5-1/2"; with drops - 4-1/4". This feature cannot be combined with a thoracic or lumbar drop. The AirHA90C comes standard with an air-controlled breakaway and pressure foot pedal.

## ▶ HEADPIECES

Add functionality to your headpiece by combining several headpiece options including Tilt, Dual-Drop, Raised, Flexion, Axial-Motion, Lateral and Translation (CBP). The AirHA90C comes standard with a Tilt and an Air Dual-Drop combination headpiece.

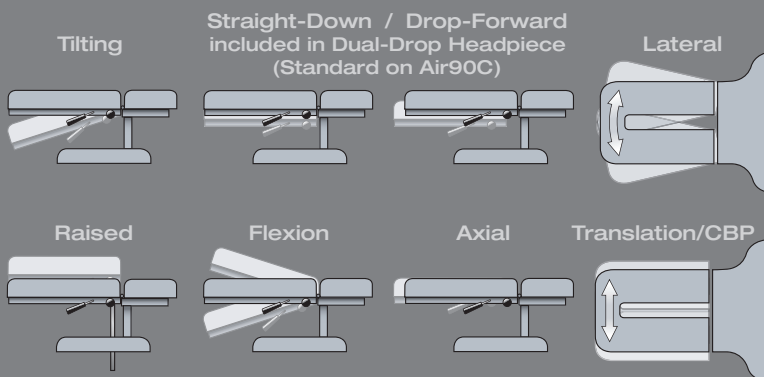


## ▶ MANUAL CERVICAL FLEXION




The Flex-Headpiece comes with the raised headpiece, occipital strap and removable tiller bar.

The headpiece unlocks and a gas-spring gently compresses and recoils as you control the flexion, lateral and long-axis movement.

All of the options below can be added to the 90C already included on the AirHA90C as standard features. These are indicated by the blue "Air" symbol.  See AirHA90C page for a complete list.

## Thoracic/Lumbar Options

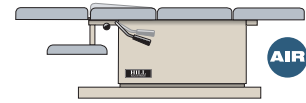


### Thoracic Drop

The HA90C manual drop is activated with an easy-lift handle. The AirHA90C is activated with Foot Strips.




### Thoracic Release

Hinges and locks at any angle down to 5 1/2 inches.

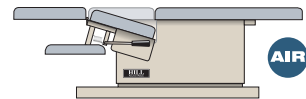

### Thoracic Breakaway

Allows for a thoracic recoil movement and features adjustable tension. The AirHA90C comes standard with an air-controlled breakaway and pressure pedal.

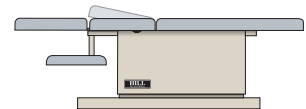

### Raised Thoracic

Tilts and locks up to 6 inches. (not available with spring breakaway or AirHA90C)



### Lumbar Drop

The HA90C manual drop is activated with an easy-lift handle. The AirHA90C is activated with Foot Strips.

OC. Many are  
features.  
AIR  
list of features.


# HA90C

More Options

Ask about our flexible lease-to-own options available for 1, 2 and 3 year financing. Your HA90C may also be eligible for a tax credit totaling nearly 1/2 the purchase price.


Visit [HillLabs.com](http://HillLabs.com) for details.

## Pelvic Options

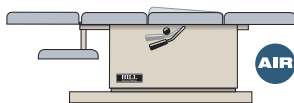


### Pelvic Drop

The HA90C manual drop is activated with an easy-lift handle. The AirHA90C is activated with Foot Strips.

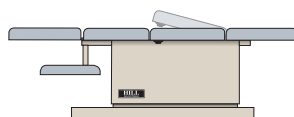

### Raised Pelvic

Tilts and locks up to 5 inches.

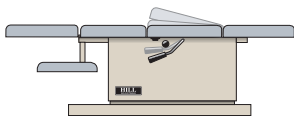

### Raised Pelvic/Drop Combo

Combines the Raised Pelvic feature with a Pelvic Drop.

*The standard pelvic cushion is not notched; optional notch is available. All HA90C drops are activated with an easy-lift handle on both sides of the table. All AirHA90C drops are selected using the control panel and activated using the Foot Strips. An opposite-side control panel is also available.*

## Other Options



### Drop-Away Foot

9" foot section folds down for treating shorter patients. Standard on AirHA90C.




### Slide-Out Ankle Support

Adjusts up to 8" in height and extends out 12" from the end of the table.




### Adjustable Armrests

Elevating armrests for supine shoulder support - often used for massage.

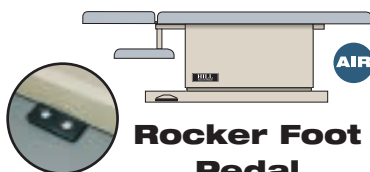

### Rocker Foot Pedal

The Rocker Foot Pedal attaches to the base on both sides of the table to control the height. Standard on the AirHA90C.

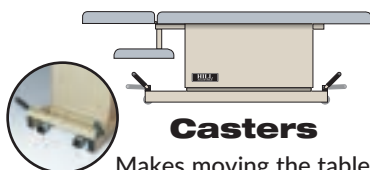

### Casters


Makes moving the table easy. Engaged with a foot lever.

## Special Configurations



### Extended Face Cut-Out

Extended face cut-out for treating children and shorter patients. Not available on AirHA90C. Other custom face cut-outs available.

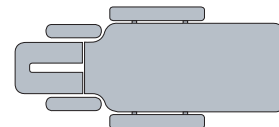


### Removable Width Extensions

3" x 21". Adds 6" to the width of the table. Used for massage therapy; often combined with adjustable armrests.



### Activator\*-Style Top

36" face cut-out and wedged foot section. 27" width, optional 30". Not available on AirHA90C.

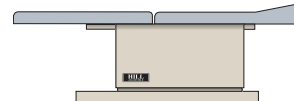

### Activator\*-Style Top

Side view showing the wedged foot section - Model HA90AM.

*The Activator\*-Style top is also available as a Hi-Lo model with soft start-up and landing, plus automatic fold-down standing board.*



# AIRHA90C

Air-Drop Table

**Crisp, smooth and impressively fast.**

Equipped with the Hill Air-Drop System\*, the AirHA90C allows doctors to select one or up to four drops simultaneously from the convenient, side-mounted, control panel. Cushions are activated instantly by tapping one of the foot strips.



*Shown in Cinder Upholstery with Black Base and added options.*

*\* The Hill Air-Drop System, with air-powered features, is only available on the AirHA90C.*

*Dual Power Foot Strips Reset the Drops with a tap!*

*Select one or up to four drops simultaneously from the Control Panel. Optional dual-sided control box available*

## ▼ STANDARD FEATURES

- Height - 21 1/2" to 31" (higher starting heights available)
- Rocker Foot Pedal and Two Power Foot Strips
- Width 24", Length 6'3"
- Armrests, 13" Face cut-out, Paper Roll, Cutter
- Tilting Headpiece, 30° positive and negative tilt
- Air Cervical Dual-Drop Headpiece
- Air-Powered Thoracic Spring Breakaway & Release
- Air Thoracic, Lumbar and Pelvic Drops
- Drop-Away Foot Section
- Choice of 20 standard colors (or upgrade to more durable Finesse upholstery)
- Beige, Gray or Black Base, enclosed for safety
- One-year warranty on all parts and labor



# SPECIFICATIONS


- Powered by an electric actuator. 110 volts, 60 Hz, 5.0 amps.; AirHA90C - 2.5 amps. 220 Hz - 50 cycle available at additional charge
- Lifting capacity - 450 lbs.
- Travel time from lowest to highest height - 15 seconds
- Choice of colors for upholstered top (see below)
- Beige, gray or black telescoping base skirting constructed of fire-rated acrylic thermo-plastic
- Powder-coated internal and external parts
- Medium density foam top - 2-1/2"; firm and soft density available
- Shipping weight - approximately 250 - 300 lbs.
- Standard width - 24"; optional widths, 22", 27", or 30"
- Range of height  
HA90C from 22" to 35"  
- with thoracic release or breakaway option, 24" to 35"  
- optional low height of 18" to 32-1/2" available  
AirHA90C from 21-1/2" to 31"  
- 23-1/2" to 33" available. Special heights upon request
- Table Top length  
HA90C - 5'10"; with drop-away foot becomes 5'1"  
AirHA90C - 6'1"; with drop-away foot becomes 5'4"  
Slide-out foot section: removed; 5'1", inserted 5'4", extended 6'4"
- Warranty - one year on all parts and labor

\*The option shown for use with the ACTIVATOR technique is designed to serve the doctor who uses the ACTIVATOR® instrument or applies the "ACTIVATOR METHODS™" technique. ACTIVATOR and ACTIVATOR METHODS are trademarks of Activator Methods, Inc. and are used for identification purposes only. Equipment made or distributed by Hill Laboratories is not made or approved by Activator Methods, Inc.

## VERY AFFORDABLE

Finance your Hill Laboratories table through our convenient and affordable lease/purchase plan. Visit HillLabs.com to calculate leasing options or call to speak personally with one of our representatives who will be happy to discuss payment options.

## AVAILABLE COLORS

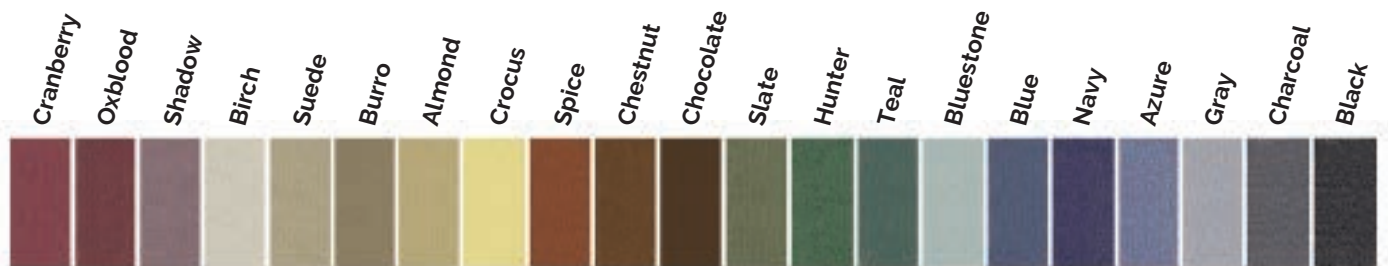

**Below - Finesse® Colors. Upgraded material; foam-backed and puncture-resistant.**



### Base Colors



Upholster your table in durable Boltaflex® with Prefixx™, scientifically formulated to give superior protection from staining and abrasions.

Finesse and custom colors are available at additional cost. Shades may vary from actual material.

Color choices may change without notice.



**Hill Laboratories**  
COMPANY

3 N. Bacton Hill Road, Frazer, PA 19355 • Toll-Free 1-877-445-5020  
610-644-2867 • fax 610-647-6297 • www.HillLabs.com