

Hill HA90P Podiatry Chair

**Beautiful. Versatile. Durable.
Affordable.**



*Shown with
Azure upholstery
and beige base*

The Hill Adjustable Podiatry Model is designed with several features usually associated with far more expensive podiatry chairs. Its quality design and craftsmanship draws from over 55 years of Hill Laboratories' manufacturing experience.

Powered by electric actuators, the HA90P has three distinct adjustable motions which position the chair smoothly, quietly, and precisely.

Adjustable height provides a comfortable 23" to 36" working range. A recline feature tilts the entire chair to any desired angle. The electric tilt back may be completely lowered to convert the HA90P to flat table.

All chair positions are controlled by convenient foot pedals. The HA90P is also available with an 'auto-return' option which automatically repositions the chair for the next patient.

Special features include a tilting foot section which can extend for taller patients, a slide-out stainless steel pan, an adjustable headrest for patient comfort and fold-down armrests. Optional features include a flexible-arm lamp and adjustable curtain (see other side).

The contour seat and back add style and comfort. The attractive vinyl upholstery is available in several colors and designed for years of continual use.

- See other side for more important information -

HA90P Features



Adjustable Armrests

Fold back out of the way when seating the patient




Adjustable Headrest

Adjusts in height 4"
Adjusts in length 12"



Adjustable Foot Section

Locks at any angle down to 85°
Extends 9.5" for taller patients.
Optional electric foot

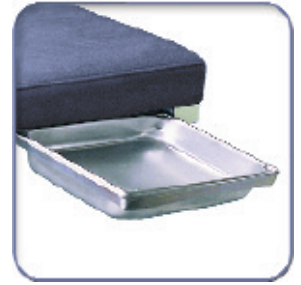

Deep Stainless Steel Pan

Removable pan slides out from under foot section



Electric Foot Pedals


Adjust the height, tilt and back



Optional Feature

Flexible Arm Lamp

Removable; comes standard with 110 outlet in base.

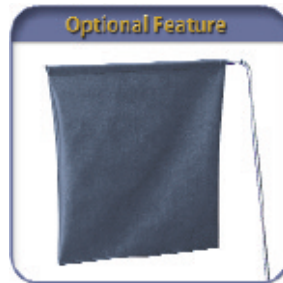

Optional Feature

Curtain

Adjusts in height




Optional Feature

Auto-Return Switch

Automatically repositions the chair for next patient

90P Positions

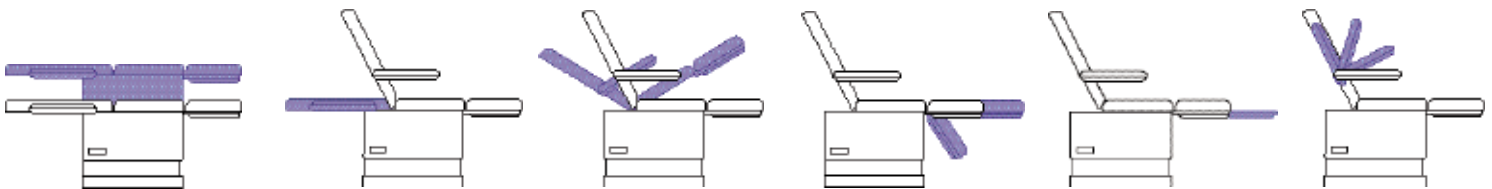

Specifications

- Elevates from low to high in 16 seconds.
- Acrylic thermo-plastic base in beige or gray.
- Durable vinyl upholstery in choice of 20 colors.
- 2 1/2" thick medium-density foam for comfort and shape retention.
- Lifting capability - 500 lbs.
- Electric requirement - 110 volts, 60 Hz, 2.7 amps.
- One-year warranty.
- Special heights available upon request.
- Shortest length when reclined - 6 feet.


Hill Laboratories
COMPANY